


DEXMED inc.

2008 Rt. 37 East, Suite 13
Toms River, NJ 08753
(732) 831-0507 (800) 833-9633
FAX (732) 831-1607

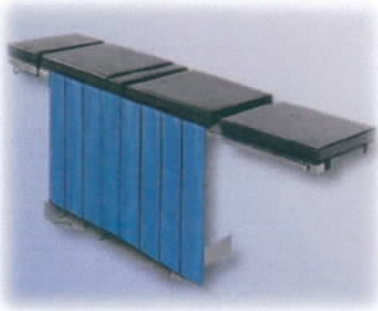

ONE SIDED VELCRO TABLE SHIELD

Available in either a single shield or a shield made with flaps. Both are custom made according to your specifications. The reverse side has a 2" strip of soft velcro which will attach to a supplied 2" strip adhesive hook. A 2" no lead strip is sewn between the velcro and the lead to allow for better draping over the table edge.

ITEM #CSUL100

SADDLE SHIELD

The saddle shield is a custom shield made for your table. The fabric that covers the table does not have any lead in it. This fabric is attached to a 2" non-lead flap to allow for better draping over the table edge. The flaps are sewn to the 2" non-lead strip and are as long as you wish. Usually we make the flaps 6" wide and they overlap 1" for added protection.

ITEM #CSSUL



SINGLE PIVOT SCATTER SHIELD

This shield is attached to a table railing with a clark socket. The socket will fit any railing up to 1" wide. A stainless pole is inserted into the shield and the short end is inserted into the socket. By releasing the handle the shield can be pivoted 180 degrees and placed wherever the operator wishes. Both a solid shield and one with flaps is available. Length is optional but width is always 36"

Pass through #805536



Solid Shield #800036